

X-LRQ-E Series Datasheet



- * 75, 150, 300, 450, 600 mm travel
- * 100 kg load capacity
- * Up to 840 mm/s speed and up to 300 N thrust
- * Built-in controller; daisy-chains with other Zaber products
- * Integrated, 500 CPR, motor mounted encoder provides slip/stall detection and recovery
- * Inline and parallel drive configurations
- * Ball screw and lead screw configurations
- * Custom versions available


Zaber's X-LRQ-E Series devices are computer-controlled, motorized linear stages with high stiffness, load, and lifetime capabilities in a compact size. Each device is available in either an inline or parallel drive configuration. They are stand-alone units requiring only a standard 24 V or 48 V power supply. The built-in motor encoder allows closed-loop operation and slip/stall recovery features. An optional indexed knob provides convenient manual control for versatile operation even without a computer.

These stages connect to the RS-232 port or USB port of any computer, and they can be daisy-chained with any other Zaber products. The daisy-chain also shares power, making it possible for multiple X-Series

products to share a single power supply. Convenient locking, 4-pin, M8 connectors on the unit allow for secure connection between units.


At only 36 mm high, these stages are excellent for applications where a low profile is required. The X-LRQ-E's innovative design allows speeds up to 840 mm/s and loads up to 100 kg. Like all of Zaber's products, the X-LRQ-E Series is designed to be 'plug and play' and very easy to set up and operate. These stages can bolt together into XY and XYZ systems. Some multi axis configurations may require additional accessories, please contact Zaber Technical Support to ensure the correct ones are selected.

Drawings

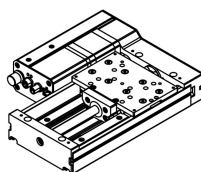

Model Number*	Travel	A	B
X-LRQ075xP-E	75.0	209.8	167.0
X-LRQ150xP-E	150.0	284.8	242.0
X-LRQ300xP-E	300.0	434.8	392.0
X-LRQ450xP-E	450.0	584.8	542.0
X-LRQ600xP-E	600.0	734.8	692.0

*See product page for complete list of available models at www.zaber.com

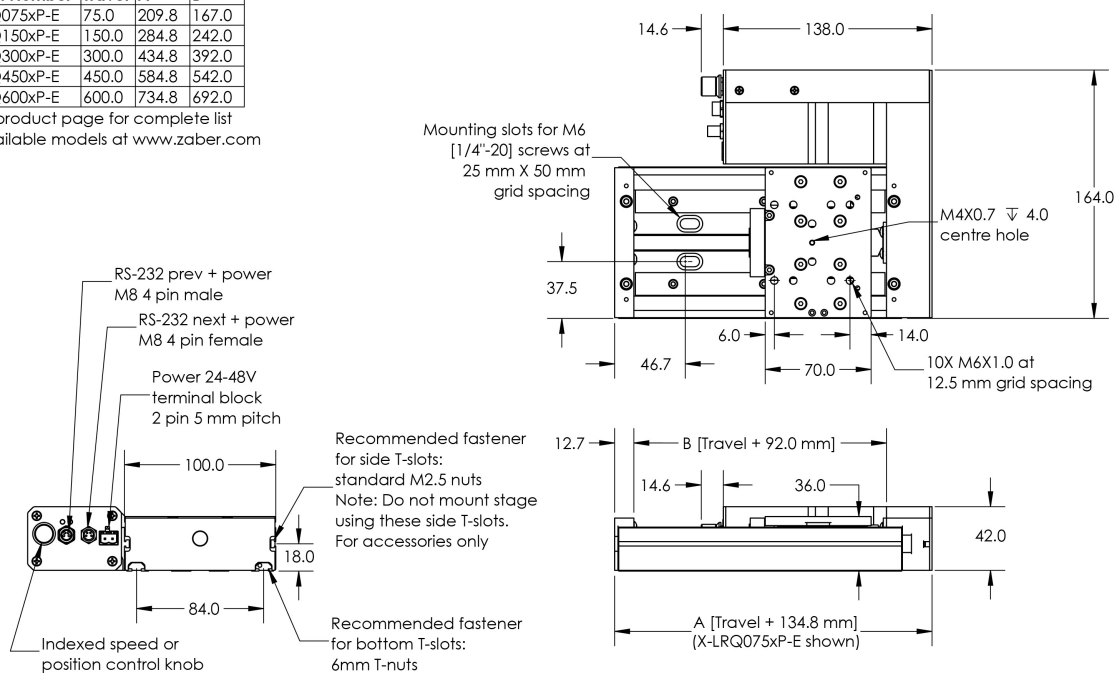

DWG 1258 R01A

Specifications

Specification	Value	Alternate Unit
Built-in Controller	Yes	
Encoder Resolution	500 CPR	2000 states/rev
Encoder Type	Rotary quadrature encoder	

Specification	Value	Alternate Unit
Communication Interface	RS-232	
Communication Protocol	Zaber ASCII (Default), Zaber Binary	
Maximum Centered Load	1000 N	224.3 lb
Maximum Cantilever Load	3000 N-cm	4248.4 oz-in
Guide Type	Recirculating Ball Linear Guide	
Stiffness in Pitch	640 N-m/°	27 µrad/N-m
Stiffness in Roll	1850 N-m/°	9 µrad/N-m
Stiffness in Yaw	665 N-m/°	26 µrad/N-m
Maximum Current Draw	1200 mA	
Power Supply	24-48 VDC	
Power Plug	2-pin screw terminal	
Motor Steps Per Rev	200	
Motor Type	Stepper (2 phase)	
Motor Rated Current	2300 mA/phase	
Inductance	2.8 mH/phase	
Default Resolution	1/64 of a step	
Data Cable Connection	Locking 4-pin M8	
Limit or Home Sensing	Magnetic home sensor	
Manual Control	Indexed knob with push switch	
Axes of Motion	1	
LED Indicators	Yes	
Mounting Interface	M6 and M3 threaded holes	
Vacuum Compatible	No	
Operating Temperature Range	0 to 50 °C	
RoHS Compliant	Yes	
CE Compliant	Yes	

Part Number	Microstep Size (Default Resolution)	Travel Range	Accuracy (unidirectional)	Repeatability
X-LRQ075AL-E01	0.09921875 µm	75 mm (2.953 ")	23 µm (0.000906 ")	< 2.5 µm (< 0.000098 ")
X-LRQ075BL-E01	0.49609375 µm	75 mm (2.953 ")	15 µm (0.000591 ")	< 2.5 µm (< 0.000098 ")

Part Number	Microstep Size (Default Resolution)	Travel Range	Accuracy (unidirectional)	Repeatability
X-LRQ075DL-E01	1.984375 µm	75 mm (2.953 ")	48 µm (0.001890 ")	< 3 µm (< 0.000118 ")
X-LRQ075HL-E01	0.1953125 µm	75 mm (2.953 ")	17.5 µm (0.000689 ")	< 5 µm (< 0.000197 ")
X-LRQ075AP-E01	0.09921875 µm	75 mm (2.953 ")	23 µm (0.000906 ")	< 2.5 µm (< 0.000098 ")
X-LRQ075BP-E01	0.49609375 µm	75 mm (2.953 ")	15 µm (0.000591 ")	< 2.5 µm (< 0.000098 ")
X-LRQ075DP-E01	1.984375 µm	75 mm (2.953 ")	48 µm (0.001890 ")	< 3 µm (< 0.000118 ")
X-LRQ075HP-E01	0.1953125 µm	75 mm (2.953 ")	27.5 µm (0.001083 ")	< 5 µm (< 0.000197 ")
X-LRQ150AL-E01	0.09921875 µm	150 mm (5.905 ")	45 µm (0.001772 ")	< 2.5 µm (< 0.000098 ")
X-LRQ150BL-E01	0.49609375 µm	150 mm (5.905 ")	25 µm (0.000984 ")	< 2.5 µm (< 0.000098 ")
X-LRQ150DL-E01	1.984375 µm	150 mm (5.905 ")	55 µm (0.002165 ")	< 3 µm (< 0.000118 ")
X-LRQ150HL-E01	0.1953125 µm	150 mm (5.905 ")	30 µm (0.001181 ")	< 5 µm (< 0.000197 ")
X-LRQ150AP-E01	0.09921875 µm	150 mm (5.905 ")	45 µm (0.001772 ")	< 2.5 µm (< 0.000098 ")
X-LRQ150BP-E01	0.49609375 µm	150 mm (5.905 ")	25 µm (0.000984 ")	< 2.5 µm (< 0.000098 ")
X-LRQ150DP-E01	1.984375 µm	150 mm (5.905 ")	55 µm (0.002165 ")	< 3 µm (< 0.000118 ")
X-LRQ150HP-E01	0.1953125 µm	150 mm (5.905 ")	40 µm (0.001575 ")	< 5 µm (< 0.000197 ")
X-LRQ300AL-E01	0.09921875 µm	300 mm (11.811 ")	90 µm (0.003543 ")	< 2.5 µm (< 0.000098 ")
X-LRQ300BL-E01	0.49609375 µm	300 mm (11.811 ")	35 µm (0.001378 ")	< 2.5 µm (< 0.000098 ")
X-LRQ300DL-E01	1.984375 µm	300 mm (11.811 ")	70 µm (0.002756 ")	< 3 µm (< 0.000118 ")
X-LRQ300HL-E01	0.1953125 µm	300 mm (11.811 ")	55 µm (0.002165 ")	< 5 µm (< 0.000197 ")
X-LRQ300AP-E01	0.09921875 µm	300 mm (11.811 ")	90 µm (0.003543 ")	< 2.5 µm (< 0.000098 ")
X-LRQ300BP-E01	0.49609375 µm	300 mm (11.811 ")	35 µm (0.001378 ")	< 2.5 µm (< 0.000098 ")

Part Number	Microstep Size (Default Resolution)	Travel Range	Accuracy (unidirectional)	Repeatability
X-LRQ300DP-E01	1.984375 µm	300 mm (11.811 ")	70 µm (0.002756 ")	< 3 µm (< 0.000118 ")
X-LRQ300HP-E01	0.1953125 µm	300 mm (11.811 ")	65 µm (0.002559 ")	< 5 µm (< 0.000197 ")
X-LRQ450AL-E01	0.09921875 µm	450 mm (17.716 ")	135 µm (0.005315 ")	< 2.5 µm (< 0.000098 ")
X-LRQ450BL-E01	0.49609375 µm	450 mm (17.716 ")	60 µm (0.002362 ")	< 2.5 µm (< 0.000098 ")
X-LRQ450DL-E01	1.984375 µm	450 mm (17.716 ")	85 µm (0.003346 ")	< 3 µm (< 0.000118 ")
X-LRQ450HL-E01	0.1953125 µm	450 mm (17.716 ")	85 µm (0.003346 ")	< 5 µm (< 0.000197 ")
X-LRQ450AP-E01	0.09921875 µm	450 mm (17.716 ")	135 µm (0.005315 ")	< 2.5 µm (< 0.000098 ")
X-LRQ450BP-E01	0.49609375 µm	450 mm (17.716 ")	60 µm (0.002362 ")	< 2.5 µm (< 0.000098 ")
X-LRQ450DP-E01	1.984375 µm	450 mm (17.716 ")	85 µm (0.003346 ")	< 3 µm (< 0.000118 ")
X-LRQ450HP-E01	0.1953125 µm	450 mm (17.716 ")	95 µm (0.003740 ")	< 5 µm (< 0.000197 ")
X-LRQ600AL-E01	0.09921875 µm	600 mm (23.622 ")	150 µm (0.005905 ")	< 2.5 µm (< 0.000098 ")
X-LRQ600BL-E01	0.49609375 µm	600 mm (23.622 ")	90 µm (0.003543 ")	< 2.5 µm (< 0.000098 ")
X-LRQ600DL-E01	1.984375 µm	600 mm (23.622 ")	100 µm (0.003937 ")	< 3 µm (< 0.000118 ")
X-LRQ600HL-E01	0.1953125 µm	600 mm (23.622 ")	115 µm (0.004528 ")	< 5 µm (< 0.000197 ")
X-LRQ600AP-E01	0.09921875 µm	600 mm (23.622 ")	150 µm (0.005905 ")	< 2.5 µm (< 0.000098 ")
X-LRQ600BP-E01	0.49609375 µm	600 mm (23.622 ")	90 µm (0.003543 ")	< 2.5 µm (< 0.000098 ")
X-LRQ600DP-E01	1.984375 µm	600 mm (23.622 ")	100 µm (0.003937 ")	< 3 µm (< 0.000118 ")
X-LRQ600HP-E01	0.1953125 µm	600 mm (23.622 ")	125 µm (0.004921 ")	< 5 µm (< 0.000197 ")
Part Number	Backlash	Maximum Speed	Minimum Speed	Speed Resolution
X-LRQ075AL-E01	< 8 µm	40 mm/s	0.000061 mm/s	0.000061 mm/s

Part Number	Backlash	Maximum Speed	Minimum Speed	Speed Resolution
	(< 0.000315 ")	(1.575 "/s)	(0.000002 "/s)	(0.000002 "/s)
X-LRQ075BL-E01	< 21 µm (< 0.000827 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ075DL-E01	< 60 µm (< 0.002362 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ075HL-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ075AP-E01	< 14 µm (< 0.000551 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ075BP-E01	< 27 µm (< 0.001063 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ075DP-E01	< 66 µm (< 0.002598 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ075HP-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ150AL-E01	< 8 µm (< 0.000315 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ150BL-E01	< 21 µm (< 0.000827 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ150DL-E01	< 60 µm (< 0.002362 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ150HL-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ150AP-E01	< 14 µm (< 0.000551 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ150BP-E01	< 27 µm (< 0.001063 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ150DP-E01	< 66 µm (< 0.002598 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ150HP-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ300AL-E01	< 8 µm (< 0.000315 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ300BL-E01	< 21 µm (< 0.000827 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ300DL-E01	< 60 µm (< 0.002362 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ300HL-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ300AP-E01	< 14 µm (< 0.000551 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)

Part Number	Backlash	Maximum Speed	Minimum Speed	Speed Resolution
X-LRQ300BP-E01	< 27 µm (< 0.001063 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ300DP-E01	< 66 µm (< 0.002598 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ300HP-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ450AL-E01	< 8 µm (< 0.000315 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ450BL-E01	< 21 µm (< 0.000827 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ450DL-E01	< 60 µm (< 0.002362 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ450HL-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ450AP-E01	< 14 µm (< 0.000551 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ450BP-E01	< 27 µm (< 0.001063 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ450DP-E01	< 66 µm (< 0.002598 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ450HP-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ600AL-E01	< 8 µm (< 0.000315 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ600BL-E01	< 21 µm (< 0.000827 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ600DL-E01	< 60 µm (< 0.002362 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ600HL-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
X-LRQ600AP-E01	< 14 µm (< 0.000551 ")	40 mm/s (1.575 "/s)	0.000061 mm/s (0.000002 "/s)	0.000061 mm/s (0.000002 "/s)
X-LRQ600BP-E01	< 27 µm (< 0.001063 ")	205 mm/s (8.071 "/s)	0.000303 mm/s (0.000012 "/s)	0.000303 mm/s (0.000012 "/s)
X-LRQ600DP-E01	< 66 µm (< 0.002598 ")	840 mm/s (33.071 "/s)	0.001211 mm/s (0.000048 "/s)	0.001211 mm/s (0.000048 "/s)
X-LRQ600HP-E01	< 30 µm (< 0.001181 ")	100 mm/s (3.937 "/s)	0.000119 mm/s (0.000005 "/s)	0.000119 mm/s (0.000005 "/s)
Part Number	Peak Thrust	Maximum Continuous Thrust	Vertical Runout	Horizontal Runout

Part Number	Peak Thrust	Maximum Continuous Thrust	Vertical Runout	Horizontal Runout
X-LRQ075AL-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075BL-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075DL-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075HL-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075AP-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075BP-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075DP-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ075HP-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 20 µm (< 0.000787 ")	< 20 µm (< 0.000787 ")
X-LRQ150AL-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150BL-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150DL-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150HL-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150AP-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150BP-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150DP-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ150HP-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 25 µm (< 0.000984 ")	< 20 µm (< 0.000787 ")
X-LRQ300AL-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ300BL-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ300DL-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ300HL-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ300AP-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 35 µm	< 30 µm

Part Number	Peak Thrust	Maximum Continuous Thrust	Vertical Runout	Horizontal Runout
			(< 0.001378 ")	(< 0.001181 ")
X-LRQ300BP-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ300DP-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ300HP-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 35 µm (< 0.001378 ")	< 30 µm (< 0.001181 ")
X-LRQ450AL-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450BL-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450DL-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450HL-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450AP-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450BP-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450DP-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ450HP-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 45 µm (< 0.001772 ")	< 40 µm (< 0.001575 ")
X-LRQ600AL-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600BL-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600DL-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600HL-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600AP-E01	110 N (24.7 lb)	100 N (22.4 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600BP-E01	100 N (22.4 lb)	100 N (22.4 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600DP-E01	25 N (5.6 lb)	25 N (5.6 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")
X-LRQ600HP-E01	300 N (67.3 lb)	200 N (44.9 lb)	< 75 µm (< 0.002953 ")	< 60 µm (< 0.002362 ")


Part Number	Pitch	Roll	Yaw	Linear Motion Per Motor Rev
X-LRQ075AL-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	1.27 mm (0.050 ")
X-LRQ075BL-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	6.35 mm (0.250 ")
X-LRQ075DL-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	25.4 mm (1.000 ")
X-LRQ075HL-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	2.5 mm (0.098 ")
X-LRQ075AP-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	1.27 mm (0.050 ")
X-LRQ075BP-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	6.35 mm (0.250 ")
X-LRQ075DP-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	25.4 mm (1.000 ")
X-LRQ075HP-E01	0.025 ° (0.436 mrad)	0.01 ° (0.174 mrad)	0.02 ° (0.349 mrad)	2.5 mm (0.098 ")
X-LRQ150AL-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	1.27 mm (0.050 ")
X-LRQ150BL-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	6.35 mm (0.250 ")
X-LRQ150DL-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	25.4 mm (1.000 ")
X-LRQ150HL-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	2.5 mm (0.098 ")
X-LRQ150AP-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	1.27 mm (0.050 ")
X-LRQ150BP-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	6.35 mm (0.250 ")
X-LRQ150DP-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	25.4 mm (1.000 ")
X-LRQ150HP-E01	0.03 ° (0.523 mrad)	0.015 ° (0.262 mrad)	0.02 ° (0.349 mrad)	2.5 mm (0.098 ")
X-LRQ300AL-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	1.27 mm (0.050 ")
X-LRQ300BL-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	6.35 mm (0.250 ")
X-LRQ300DL-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	25.4 mm (1.000 ")
X-LRQ300HL-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	2.5 mm (0.098 ")
X-LRQ300AP-E01	0.034 °	0.015 °	0.03 °	1.27 mm (0.050 ")

Part Number	Pitch	Roll	Yaw	Linear Motion Per Motor Rev
	(0.593 mrad)	(0.262 mrad)	(0.523 mrad)	
X-LRQ300BP-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	6.35 mm (0.250 ")
X-LRQ300DP-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	25.4 mm (1.000 ")
X-LRQ300HP-E01	0.034 ° (0.593 mrad)	0.015 ° (0.262 mrad)	0.03 ° (0.523 mrad)	2.5 mm (0.098 ")
X-LRQ450AL-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	1.27 mm (0.050 ")
X-LRQ450BL-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	6.35 mm (0.250 ")
X-LRQ450DL-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	25.4 mm (1.000 ")
X-LRQ450HL-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	2.5 mm (0.098 ")
X-LRQ450AP-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	1.27 mm (0.050 ")
X-LRQ450BP-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	6.35 mm (0.250 ")
X-LRQ450DP-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	25.4 mm (1.000 ")
X-LRQ450HP-E01	0.04 ° (0.698 mrad)	0.025 ° (0.436 mrad)	0.04 ° (0.698 mrad)	2.5 mm (0.098 ")
X-LRQ600AL-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	1.27 mm (0.050 ")
X-LRQ600BL-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	6.35 mm (0.250 ")
X-LRQ600DL-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	25.4 mm (1.000 ")
X-LRQ600HL-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	2.5 mm (0.098 ")
X-LRQ600AP-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	1.27 mm (0.050 ")
X-LRQ600BP-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	6.35 mm (0.250 ")
X-LRQ600DP-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	25.4 mm (1.000 ")
X-LRQ600HP-E01	0.045 ° (0.785 mrad)	0.035 ° (0.611 mrad)	0.04 ° (0.698 mrad)	2.5 mm (0.098 ")


Part Number	Mechanical Drive System	Weight
X-LRQ075AL-E01	Precision lead screw	2.19 kg (4.828 lb)
X-LRQ075BL-E01	Precision lead screw	2.19 kg (4.828 lb)
X-LRQ075DL-E01	Precision lead screw	2.19 kg (4.828 lb)
X-LRQ075HL-E01	Precision ball screw	2.19 kg (4.828 lb)
X-LRQ075AP-E01	Precision lead screw	2.33 kg (5.137 lb)
X-LRQ075BP-E01	Precision lead screw	2.33 kg (5.137 lb)
X-LRQ075DP-E01	Precision lead screw	2.33 kg (5.137 lb)
X-LRQ075HP-E01	Precision ball screw	2.33 kg (5.137 lb)
X-LRQ150AL-E01	Precision lead screw	2.56 kg (5.644 lb)
X-LRQ150BL-E01	Precision lead screw	2.56 kg (5.644 lb)
X-LRQ150DL-E01	Precision lead screw	2.56 kg (5.644 lb)
X-LRQ150HL-E01	Precision ball screw	2.56 kg (5.644 lb)
X-LRQ150AP-E01	Precision lead screw	2.70 kg (5.952 lb)
X-LRQ150BP-E01	Precision lead screw	2.70 kg (5.952 lb)
X-LRQ150DP-E01	Precision lead screw	2.70 kg (5.952 lb)
X-LRQ150HP-E01	Precision ball screw	2.70 kg (5.952 lb)
X-LRQ300AL-E01	Precision lead screw	3.32 kg (7.319 lb)
X-LRQ300BL-E01	Precision lead screw	3.32 kg (7.319 lb)
X-LRQ300DL-E01	Precision lead screw	3.32 kg (7.319 lb)
X-LRQ300HL-E01	Precision ball screw	3.32 kg (7.319 lb)
X-LRQ300AP-E01	Precision lead screw	3.46 kg (7.628 lb)
X-LRQ300BP-E01	Precision lead screw	3.46 kg (7.628 lb)

Part Number	Mechanical Drive System	Weight
X-LRQ300DP-E01	Precision lead screw	3.46 kg (7.628 lb)
X-LRQ300HP-E01	Precision ball screw	3.46 kg (7.628 lb)
X-LRQ450AL-E01	Precision lead screw	4.07 kg (8.973 lb)
X-LRQ450BL-E01	Precision lead screw	4.07 kg (8.973 lb)
X-LRQ450DL-E01	Precision lead screw	4.07 kg (8.973 lb)
X-LRQ450HL-E01	Precision ball screw	4.07 kg (8.973 lb)
X-LRQ450AP-E01	Precision lead screw	4.21 kg (9.281 lb)
X-LRQ450BP-E01	Precision lead screw	4.21 kg (9.281 lb)
X-LRQ450DP-E01	Precision lead screw	4.21 kg (9.281 lb)
X-LRQ450HP-E01	Precision ball screw	4.21 kg (9.281 lb)
X-LRQ600AL-E01	Precision lead screw	4.78 kg (10.538 lb)
X-LRQ600BL-E01	Precision lead screw	4.78 kg (10.538 lb)
X-LRQ600DL-E01	Precision lead screw	4.78 kg (10.538 lb)
X-LRQ600HL-E01	Precision ball screw	4.78 kg (10.538 lb)
X-LRQ600AP-E01	Precision lead screw	4.92 kg (10.847 lb)
X-LRQ600BP-E01	Precision lead screw	4.92 kg (10.847 lb)
X-LRQ600DP-E01	Precision lead screw	4.92 kg (10.847 lb)
X-LRQ600HP-E01	Precision ball screw	4.92 kg (10.847 lb)


Charts

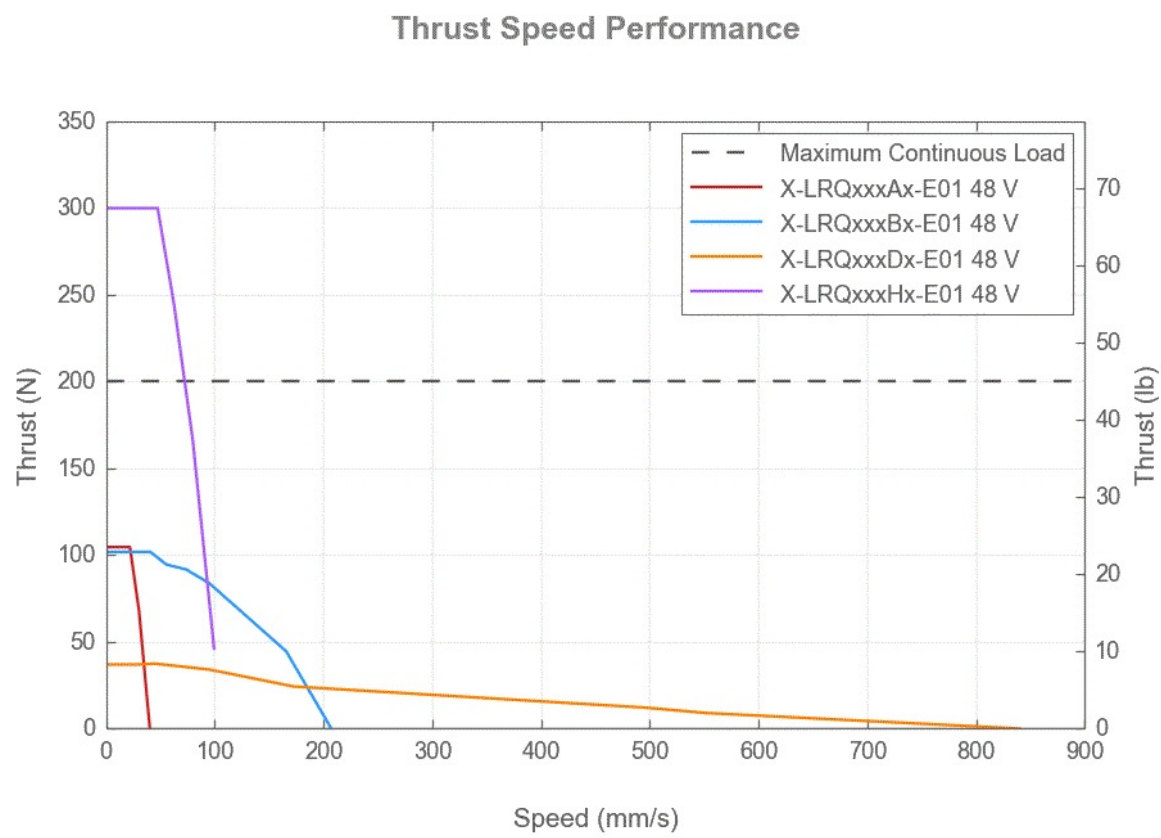


Thrust Speed Performance

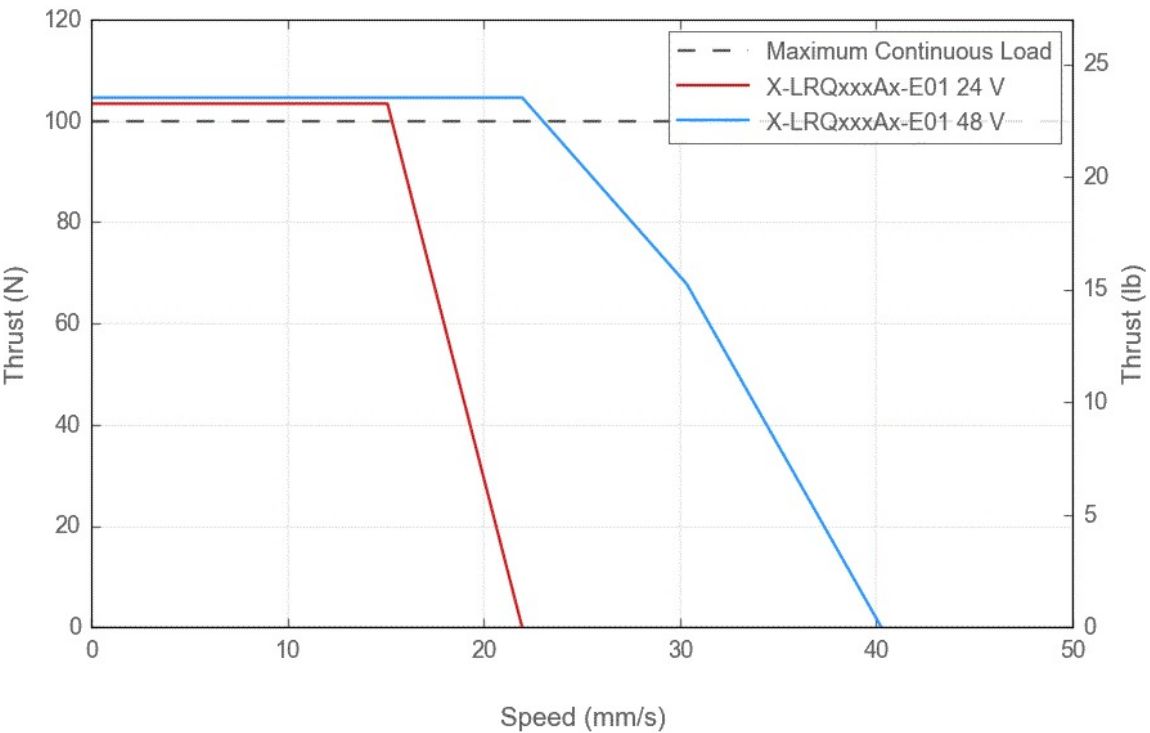


Thrust Speed Performance

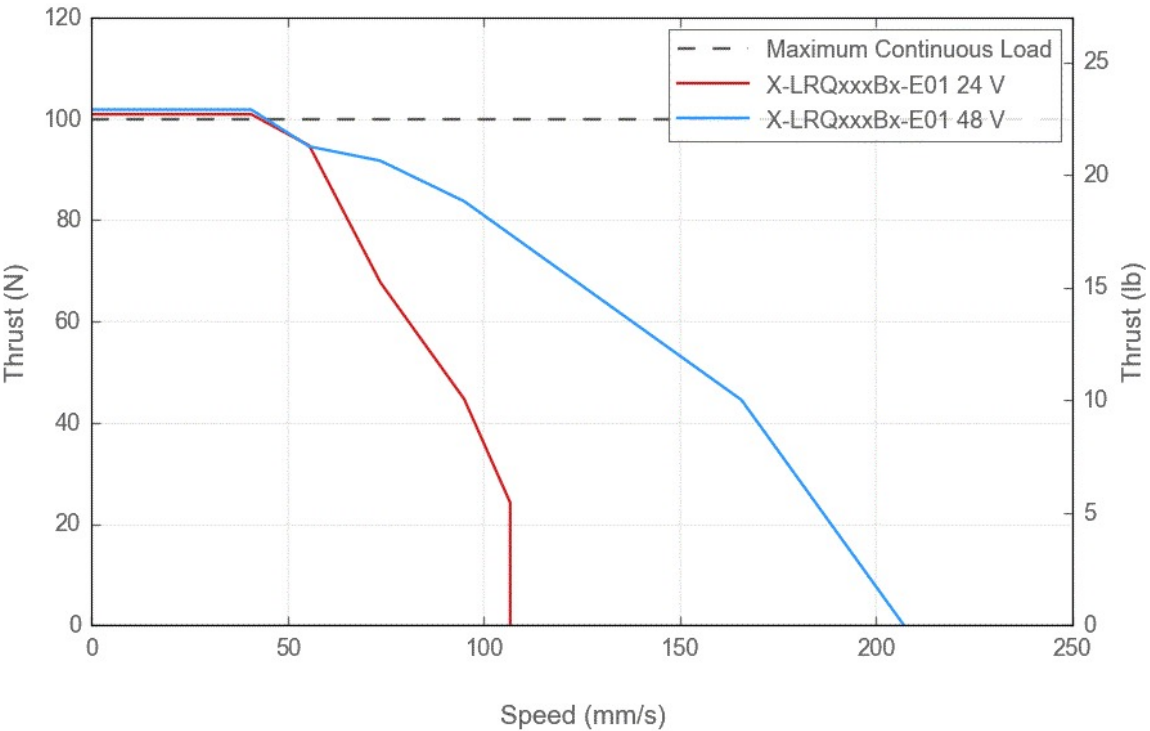

Thrust Speed Performance

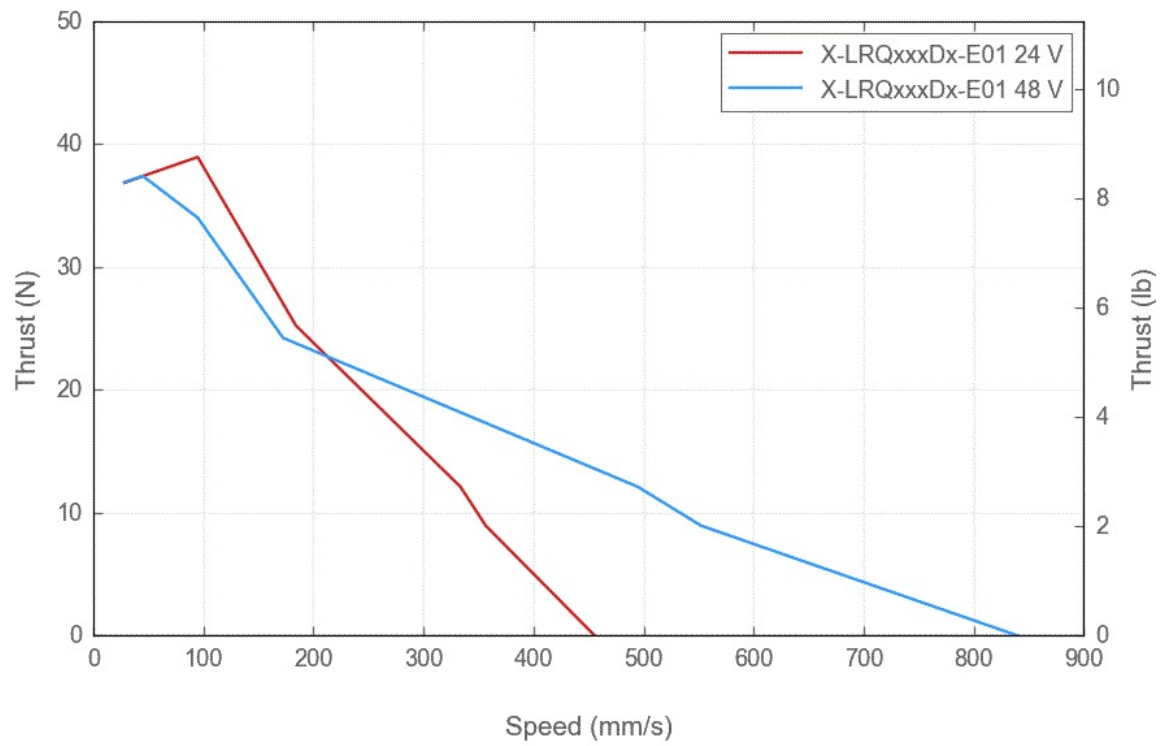

Thrust Speed Performance

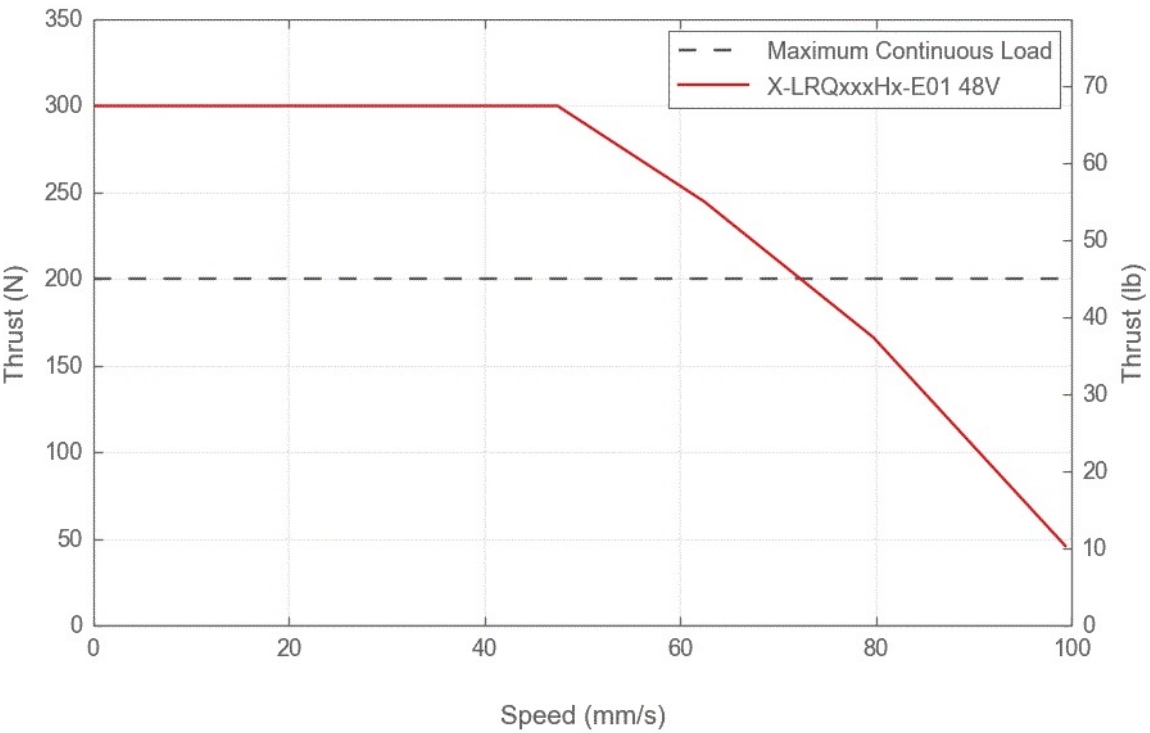

LRQ Linear Bearing Lifetime

