



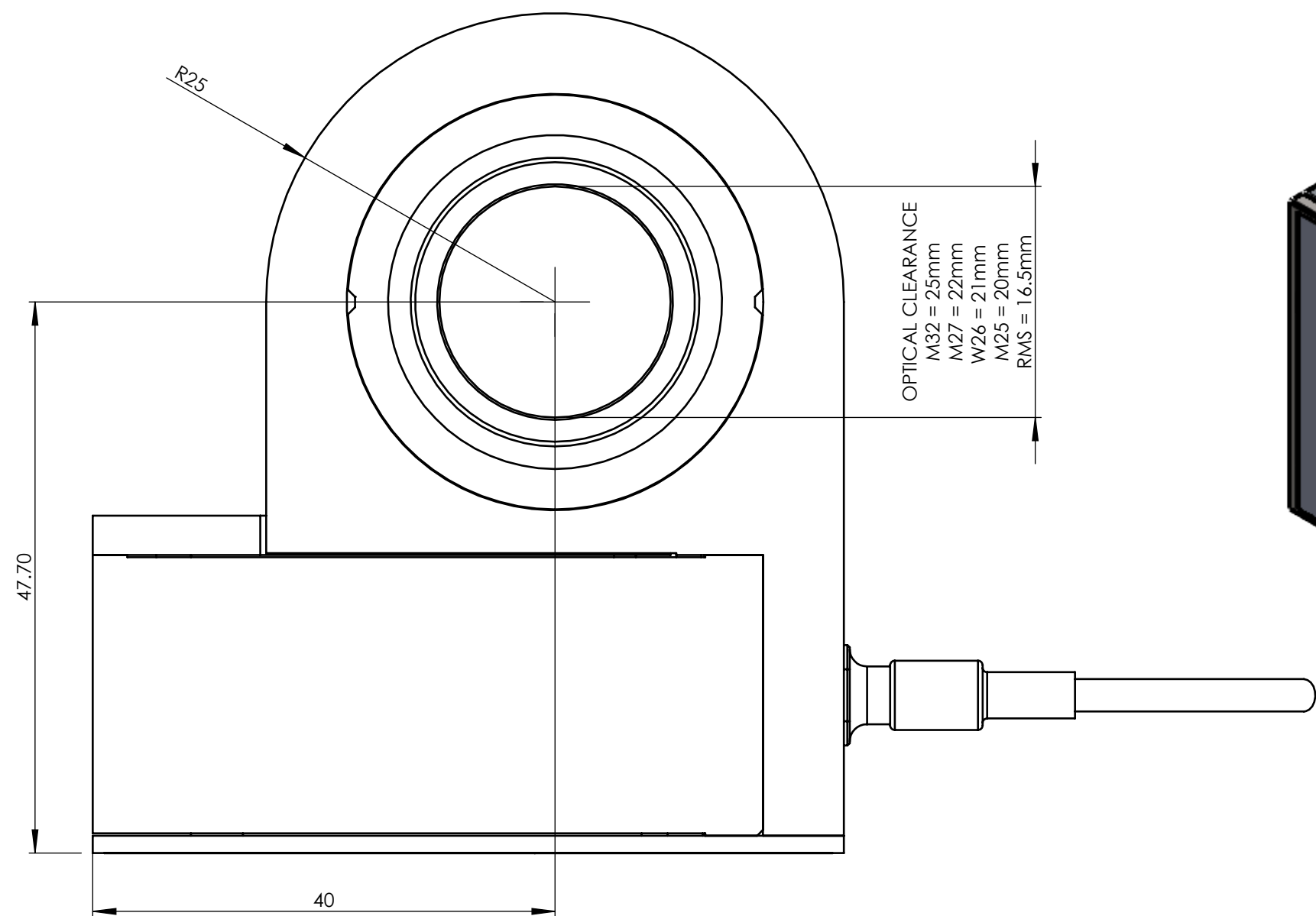


UPRIGHT ORIENTATION SHOWN

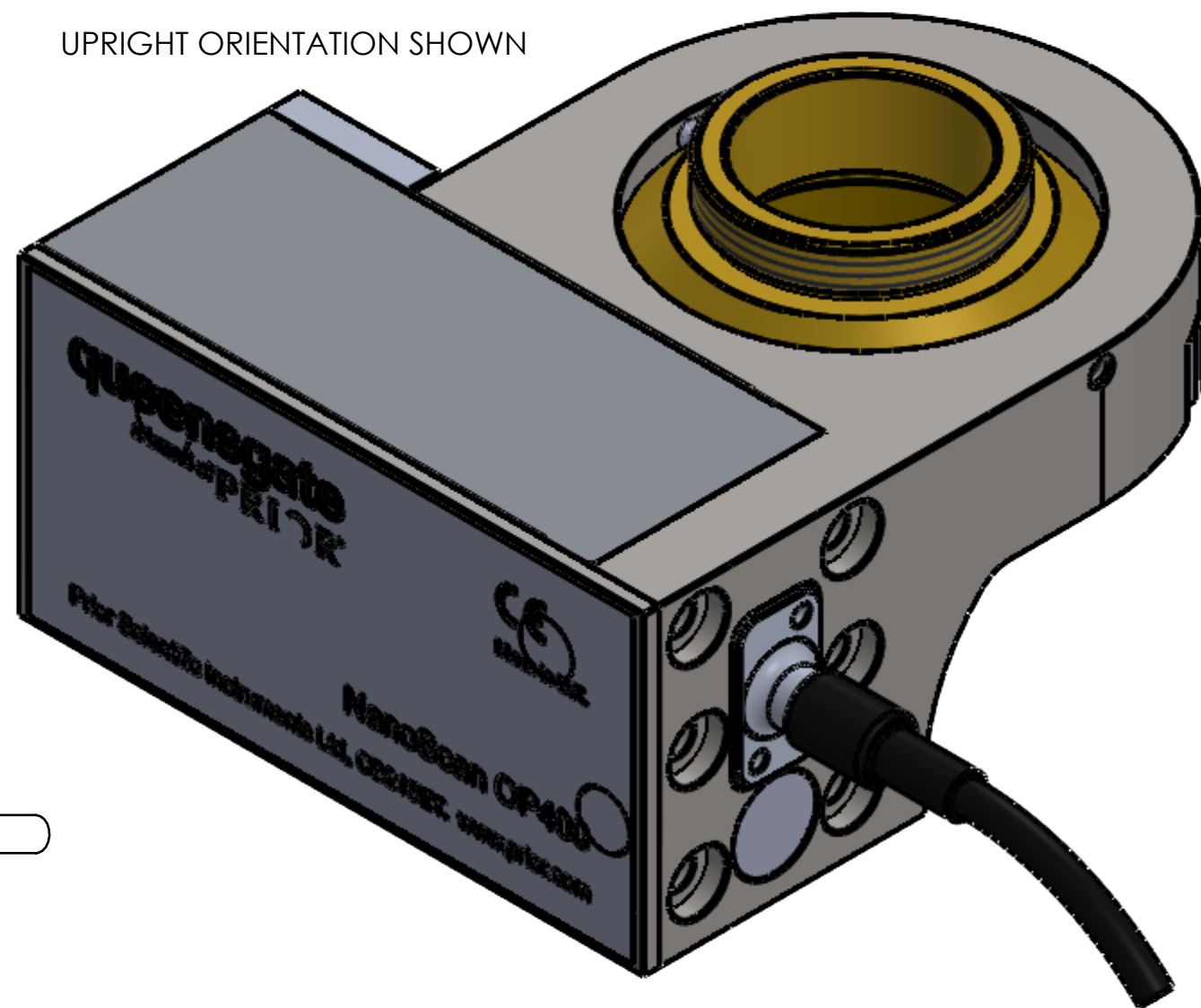

ADAPTERS AND ACCESSORIES

(SOLD SEPARATLY NOT SHOWN TO SCALE)

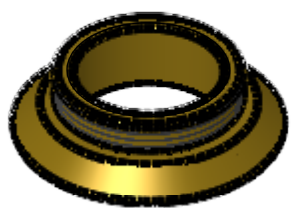


MICROSCOPE THREAD MOUNT. SELECT THREAD TO MATCH MICROSCOPE RMS, M25, W26, M27 OR M32.






OBJECTIVE THREAD ADAPTER. SELECT APPROPRIATE THREAD TO MATCH OBJECTIVE RMS, M25, W26 OR M27. NOT REQUIRED FOR M32 SYSTEMS.



PARAFOCAL SPACERS FOR MOVING OTHER OBJECTIVES IN MICROSCOPE NOSE TO THE SAME FOCAL HEIGHT RMS, M25, W26, M27 OR M32. 13mm LONG EXCEPT M32 WHICH IS 12mm.

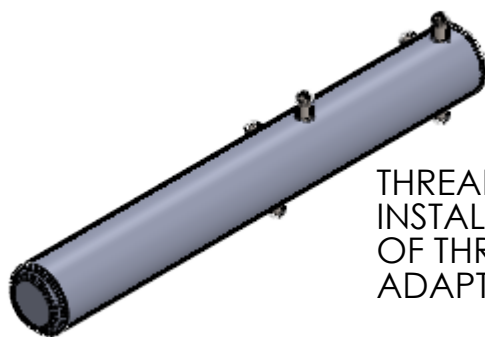

THREAD ADAPTER KEY CAN AID INSTALLATION AND REMOVAL OF THREAD ADAPTERS LIGHTLY NIP UP ADAPTERS ON INSTALLATION.

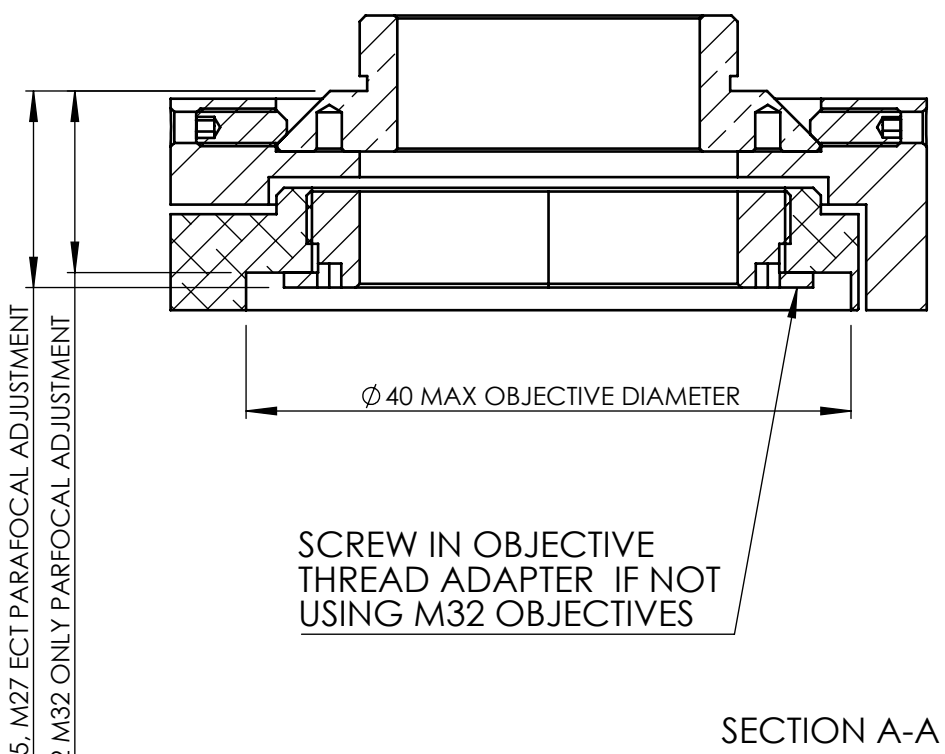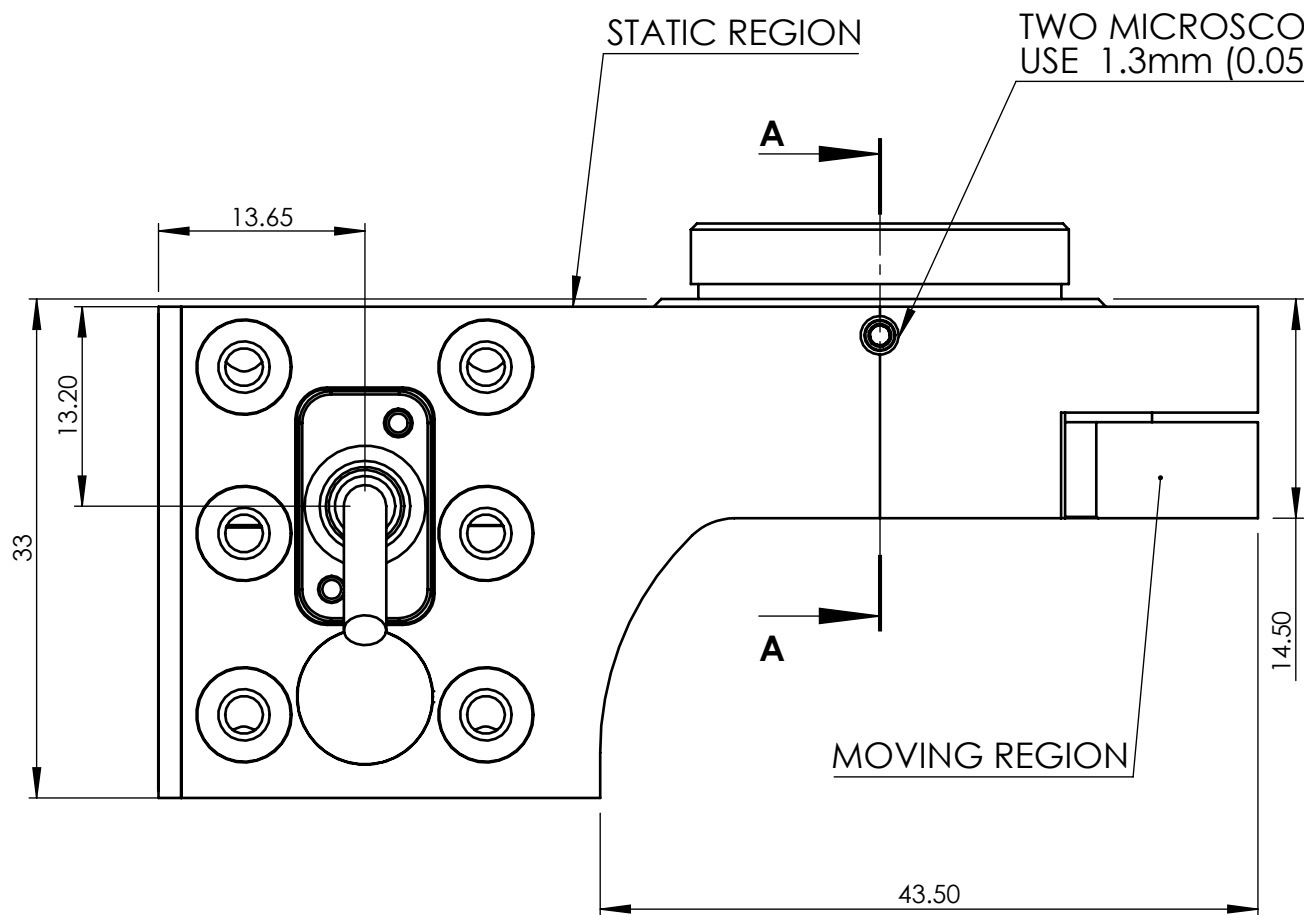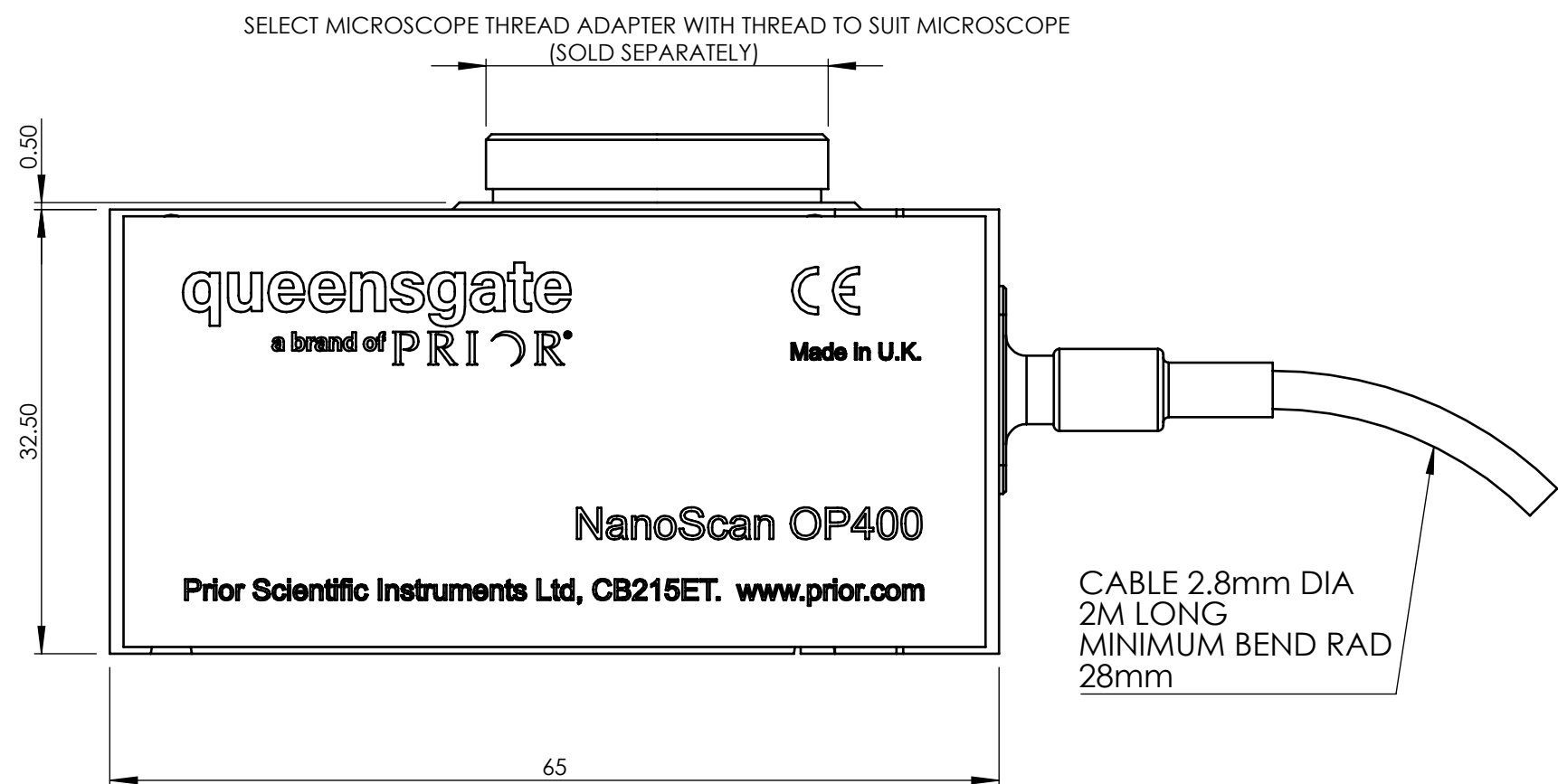

TITLE	OP-400 INSTALLATION DRAWING	APRVD.		3RD ANGLE	DO NOT SCALE IF IN DOUBT ASK
DRG. No.	QG-0494-IN-100001	DRAWN	JC		
MATL.	STAINLESS, BRASS, ALUMINIUM	SCALE	1:1	TOLERANCE UNLESS OTHERWISE STATED	
FINISH	CLEAN	DATE	09.08.19	±0.1mm	
PRIOR Scientific Instruments Ltd. Fulbourn, Cambridge.				BREAK ALL SHARP EDGES	

ECN	REV	CREATED	APP	DATE	DETAILS
1st Rel	A	JC			1st Release